
ALIGNERS

Scanner Export **and Integration** 2026

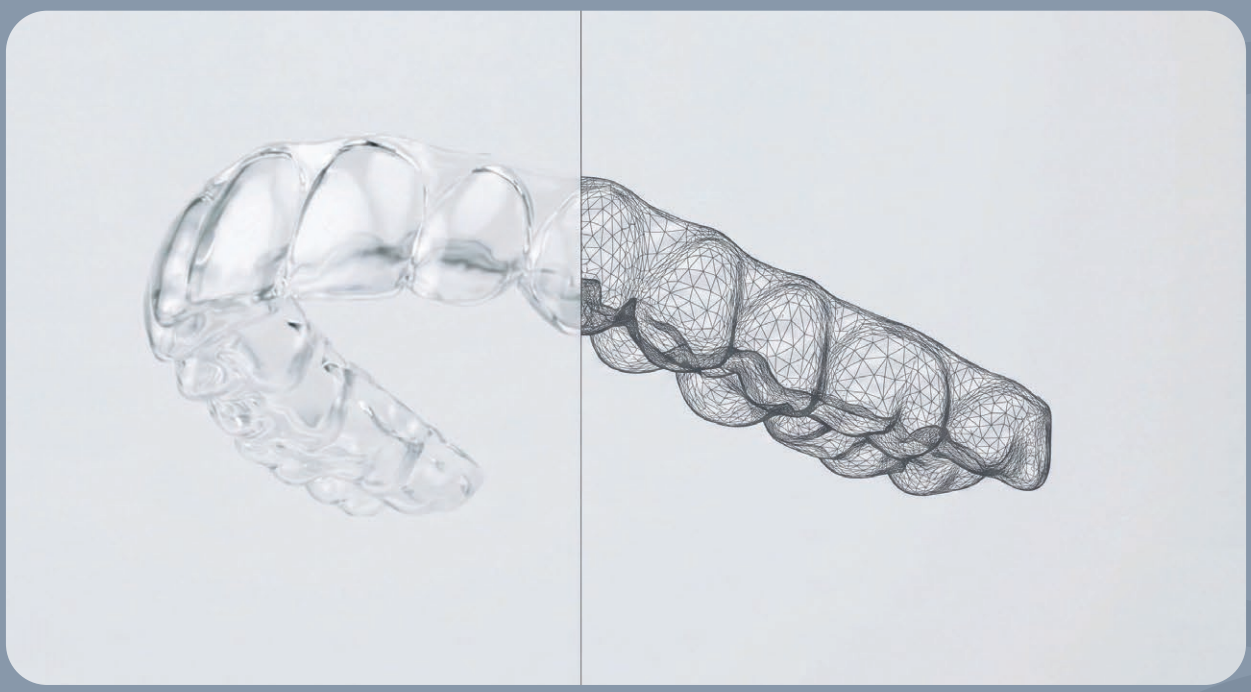

www.gloaligners.co.uk

Content

01
3Shape Integration

02
3Shape STL Export


03
Medit Integration

04
Medit STL Export

05
Itero STL Export

01
3Shape Integration

1. Head to the main menu on your doctor portal.

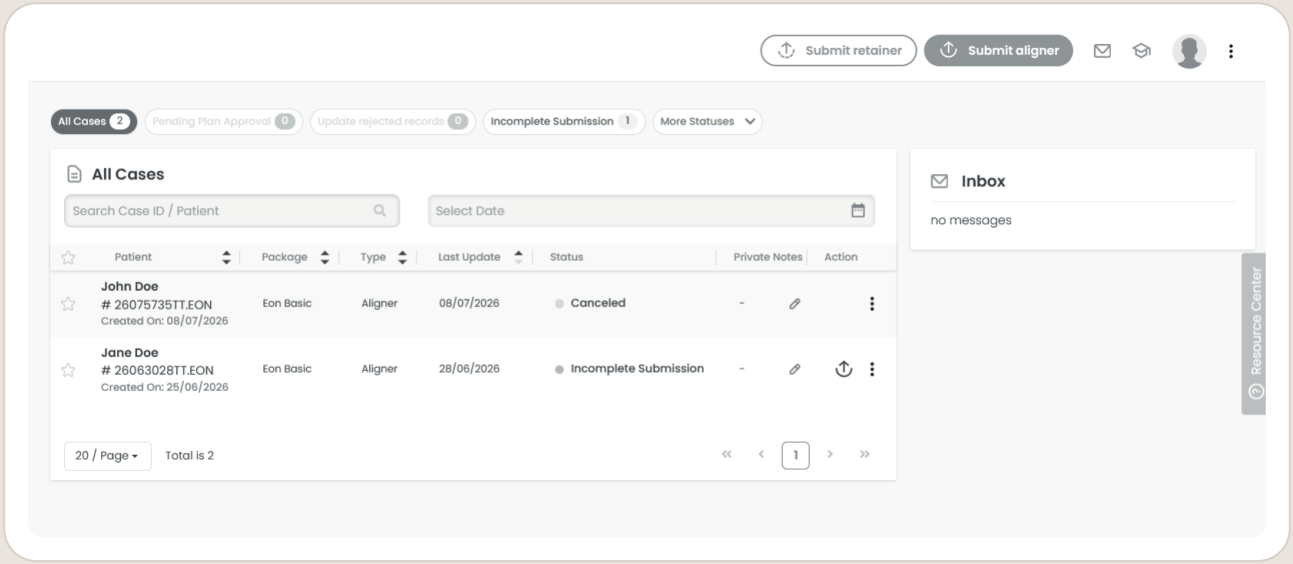


2. Go to the Integrations tab within Settings.


3. Find 3Shape scanner integration option and click Connect.




4 Enter your email address.

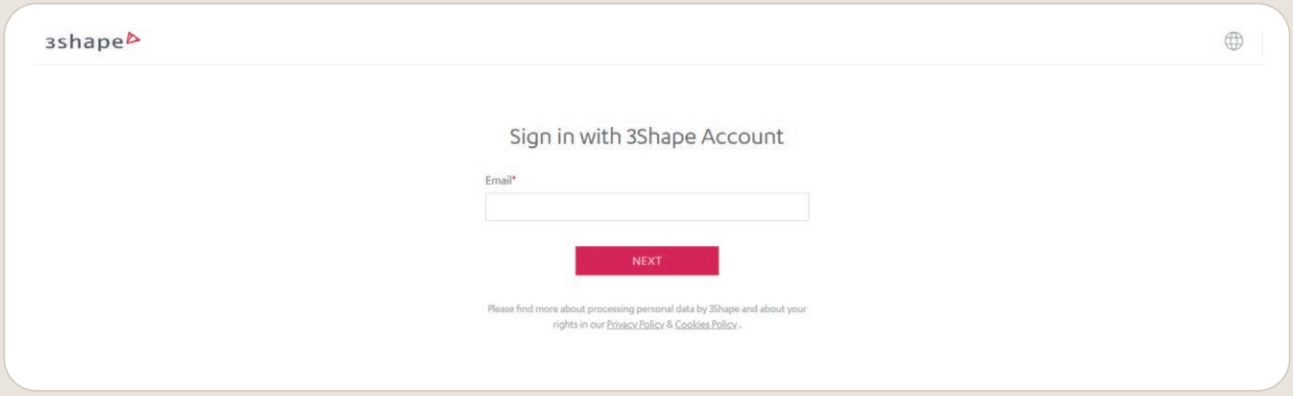

5 Enter your password.

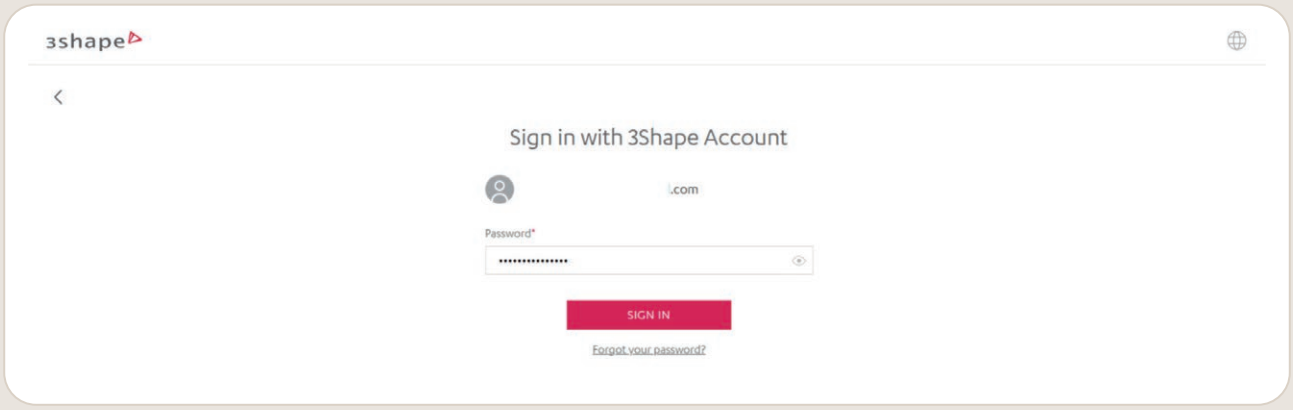

6 Your 3Shape account is now connected to your doctor portal.

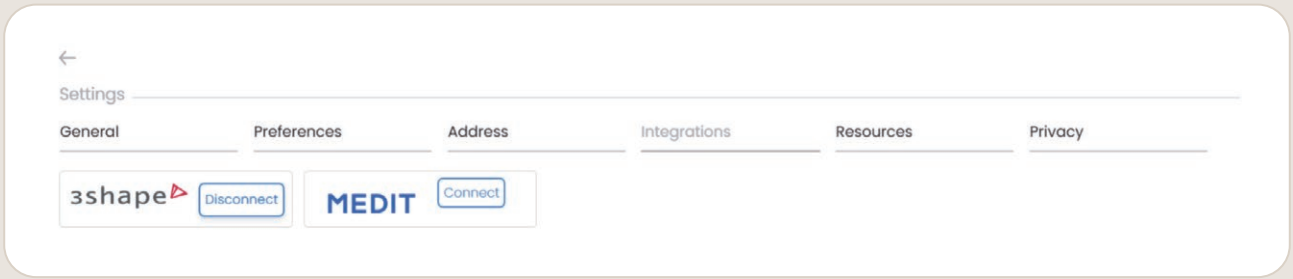

02
3Shape STL Export

Export Result Files from Case List


TRIOS scans can also be exported in STL, DCM, and PLY formats.

> **Note:** The TRIOS Scan STL Export option must be enabled on the dongle to export TRIOS scans in STL format.

If the "Export Scan" option is unavailable, contact your reseller or 3Shape Support to verify that the required license is enabled.


1. In Dental Manager, right click the TRIOS order.
2. Select "Advanced".
3. Click "Export Scan".
4. When prompted with "Save model in Y-axis or Z-axis?", select the desired orientation.
   - **Y-axis:** Standard export orientation.
   - **Z-axis:** Use this option if the destination software requires a different model orientation.


     (This option can be used if the 3rd party software for which the scans will be used corresponds another scan direction).

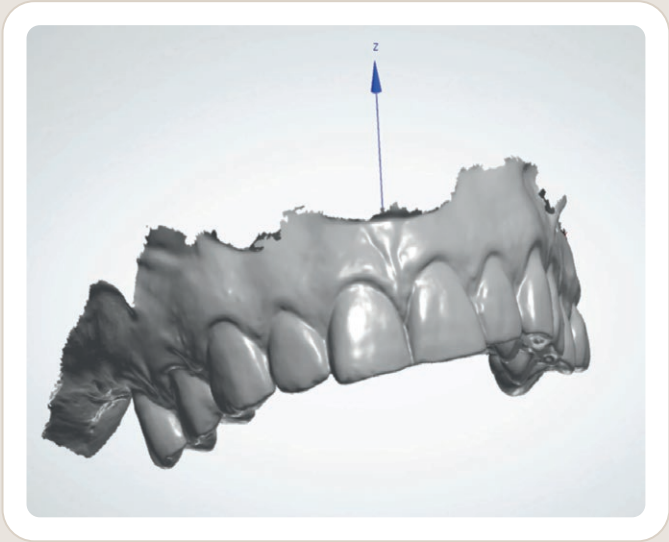

5 When prompted with "Decimate scan?", choose whether to reduce the mesh size.

- Select Yes to reduce the mesh size and specify the desired decimation percentage.
- Select No to export the scan in its original resolution.

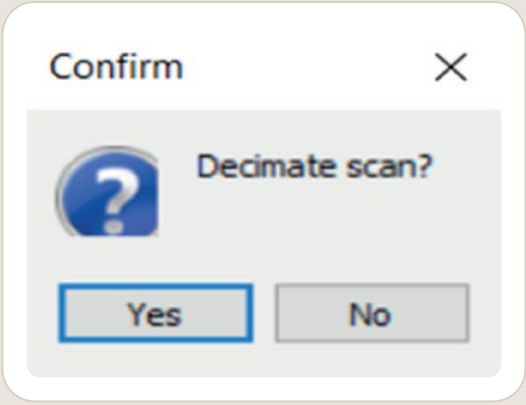


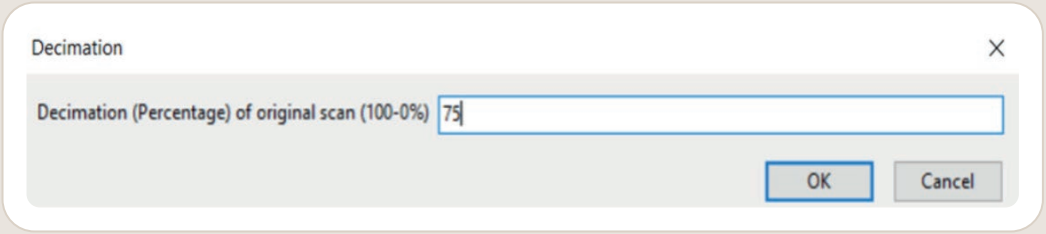

7 Choose the destination folder.


8 Click "Save" to complete the export.

Once files are exported

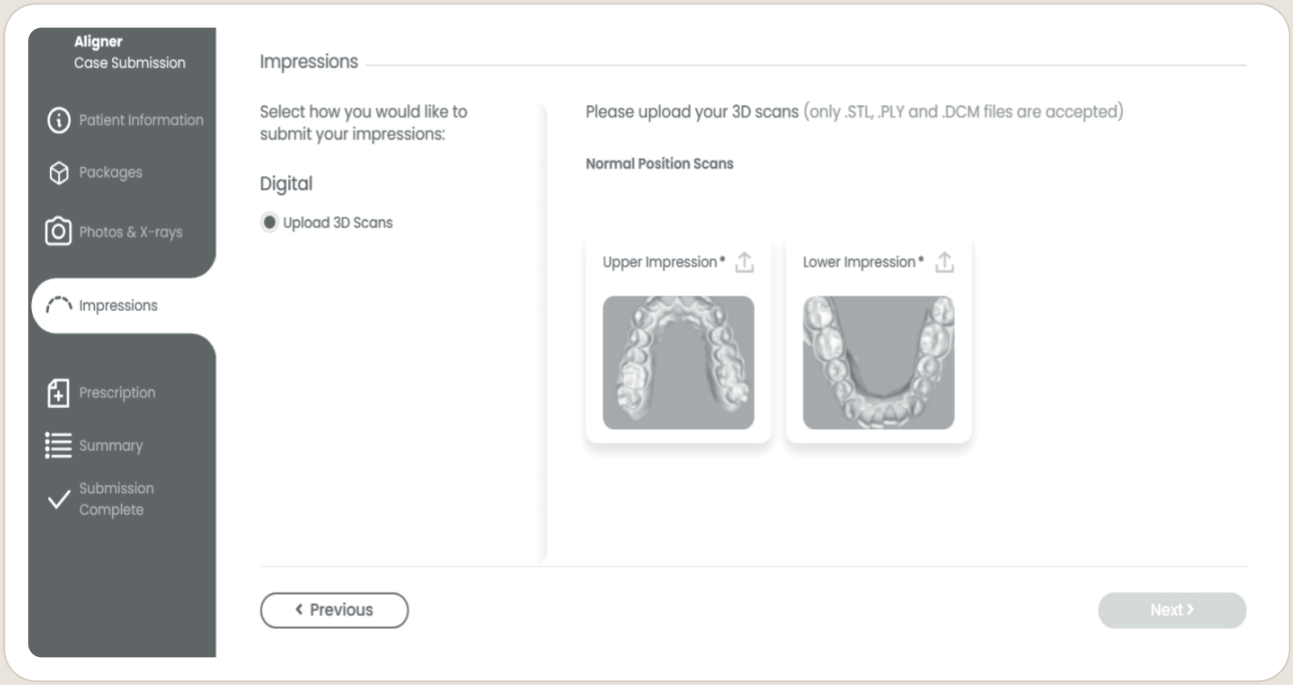


03
Medit Integration

1. Head to the main menu on your doctor portal.

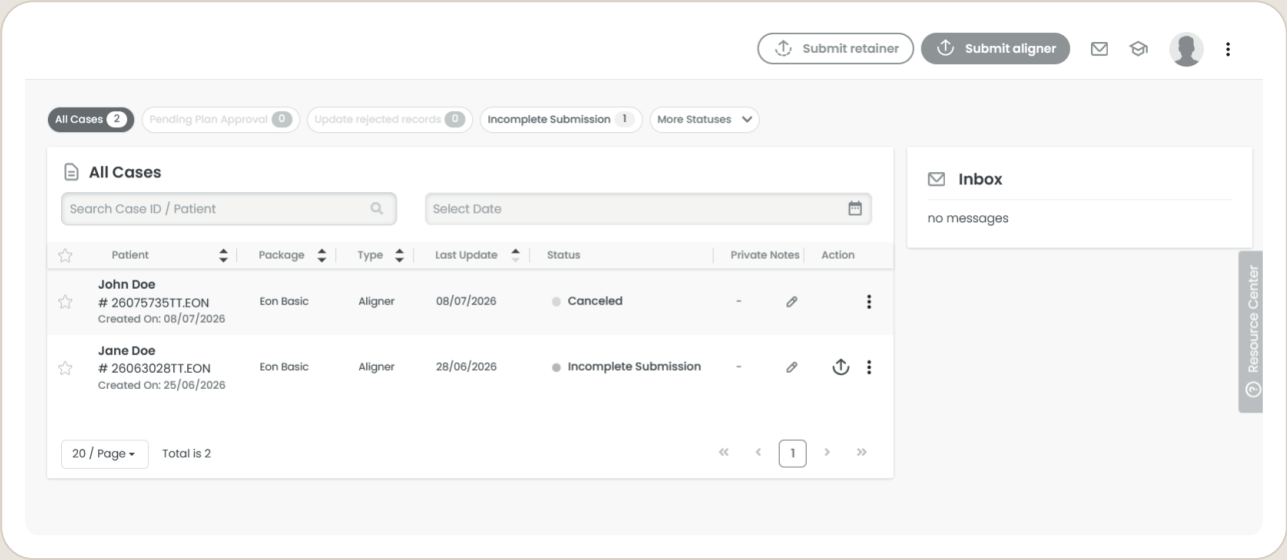

2. Go to the Integrations tab within Settings.

3. Find Medit scanner integration option and click Connect.

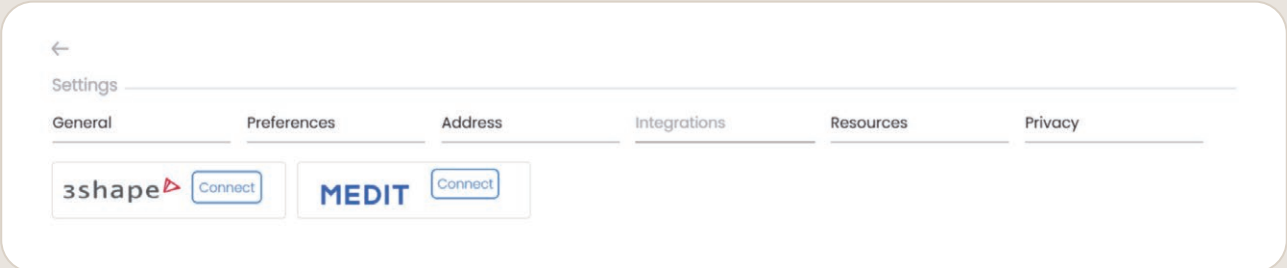

4 Enter your email and password.




5 Your Medit account is now connected to your doctor portal.




04
Medit STL Export

Export Result Files from Case List


You can export result files of a case or multiple cases simultaneously with the "Export" option above the case list in Case Box and In Box.

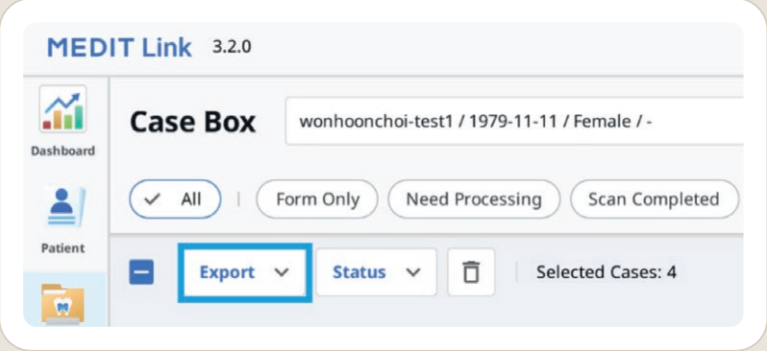


1. Go to In Box or Work Box.
2. Select all cases you want to export result files from the list.

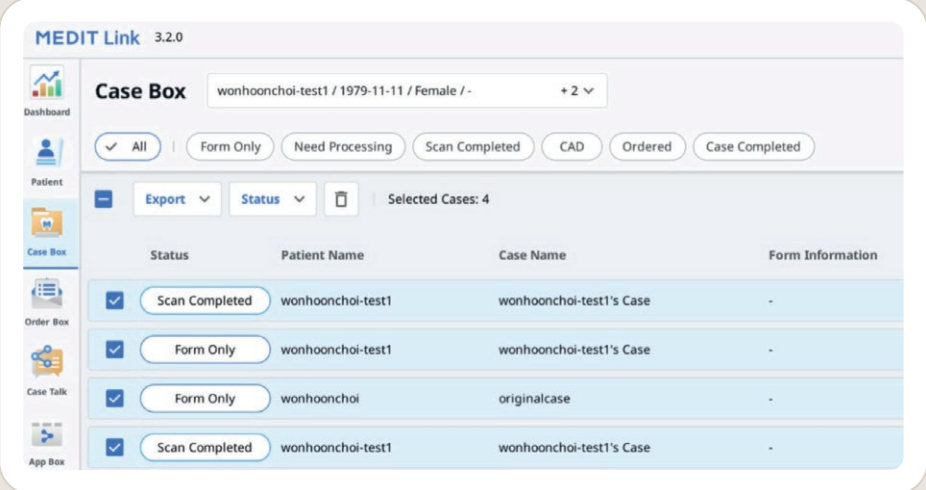

3 Click "Export" and select "Result Files".

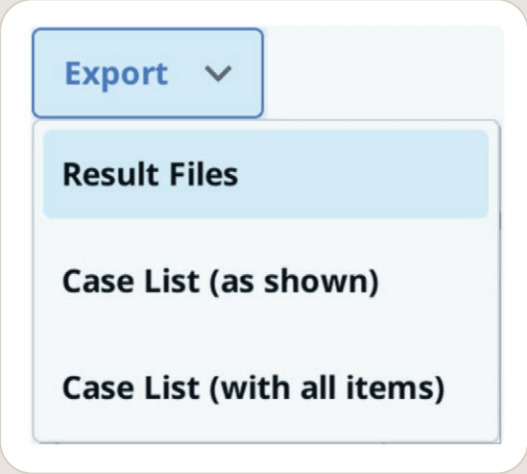


4 Set the destination folder, file type, mesh format, and image format, and select export options on the "Options" tab.

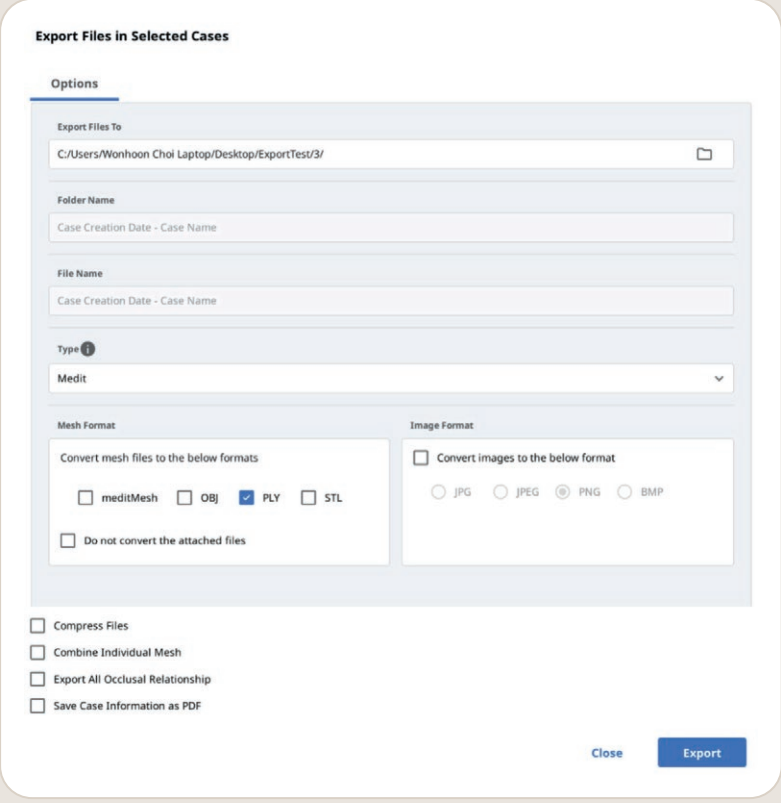


5 Click "Export".

6 Click "Yes" to open the folder where the files are saved.

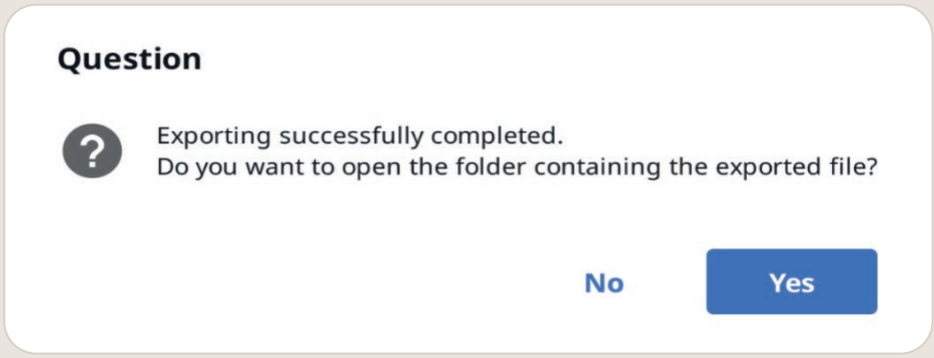

Once files are exported


After downloading your files, you can upload them by using the Upload 3D Scans from Impressions tab.

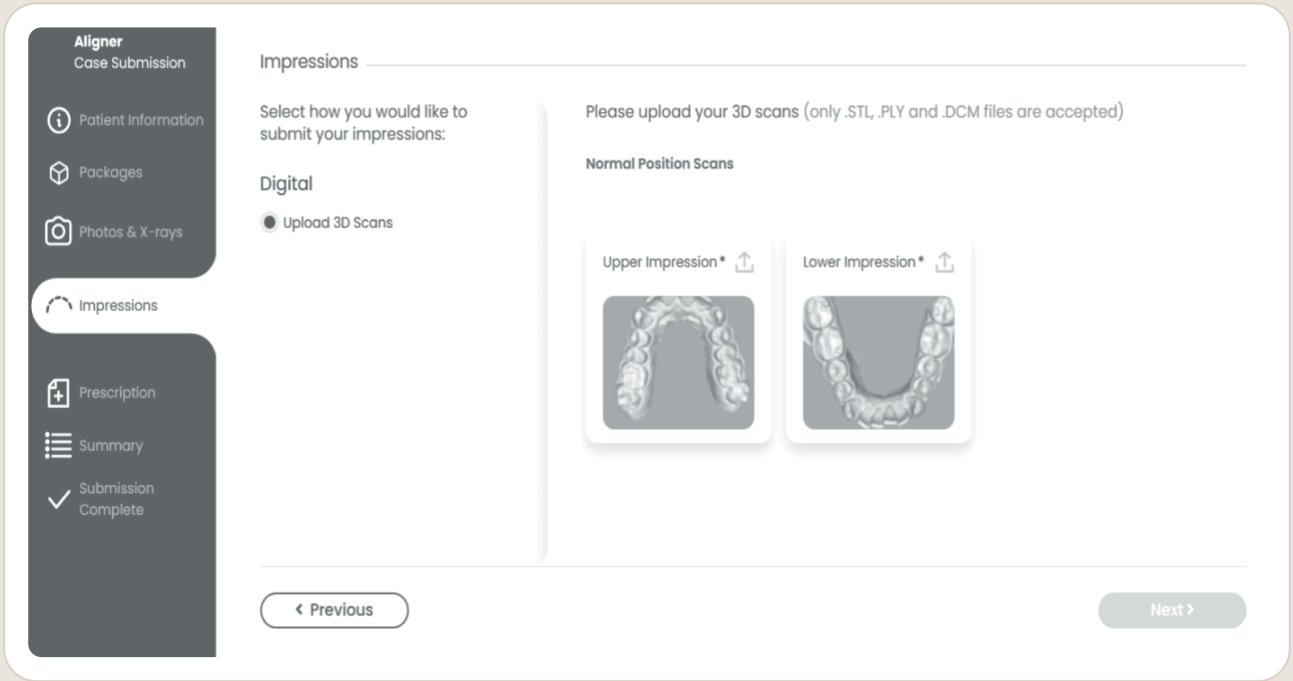

05
Itero STL Export

1. Log in to MyiTero.com™ (MIDC) on your computer.

2. Select the case that you would like to export and click 'Export'.

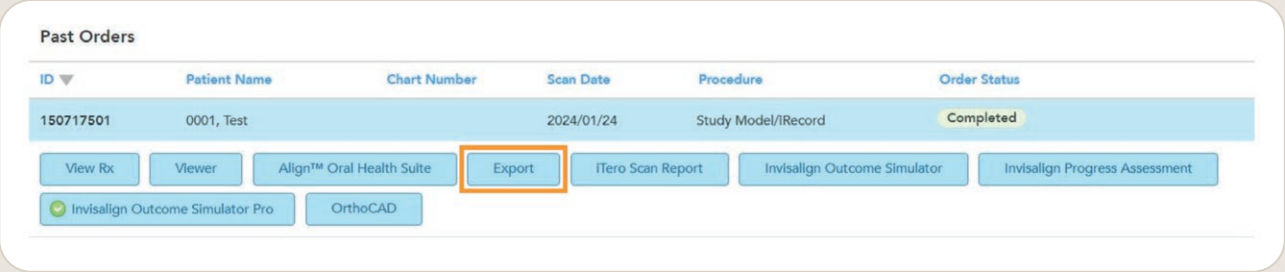


3 A. The following is the recommended export option for an Orthodontic case:

Export Type: Open shell.


Data Format: File per Arch (arches oriented in occlusion).

File Type: STL + PLY.

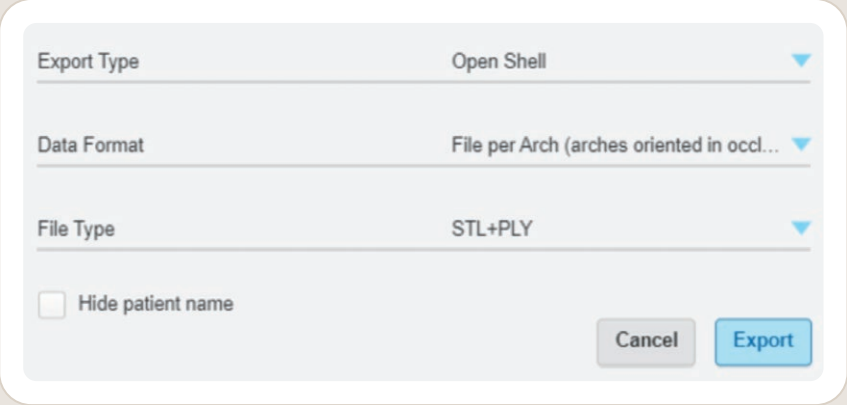


3 B. The following is the recommended export option for a restorative case:

CAD/CAM System: Implants and other exports, open model.

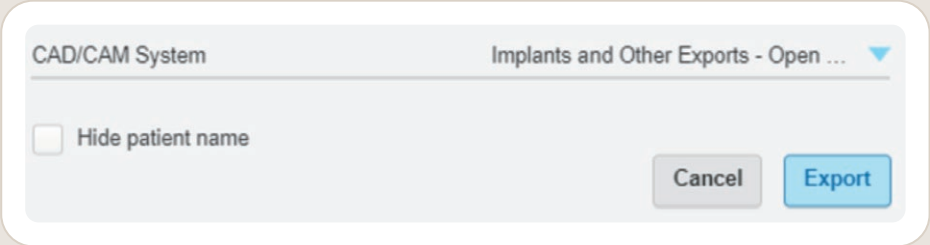

4 All exported files will be available in the computer or laptop Download folder.

Once files are exported

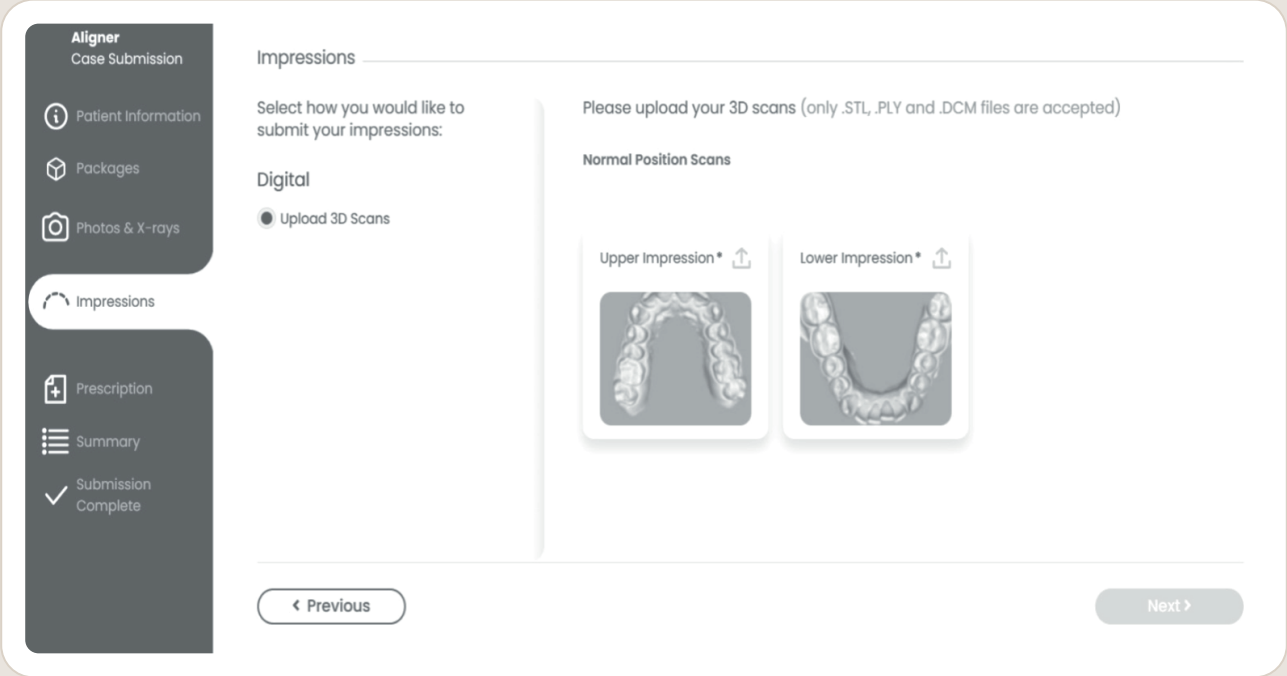

glo

ALIGNERS

Predictable outcomes begin
with predictable planning.

Glo Aligners - planned by the founders
of the Complete Aligner Programme.

gloaligners.co.uk contact@gloaligners.co.uk
Unit 22/23 Howley Park Business Village, Morley, Leeds, LS27 0BZ